# ORTOCEL Soft

ENGLISH




*Cold/hot gel pack designed to reduce pain and inflammation in post-surgical and post-traumatic treatments.*

*It is made of a gelatinous emulsion that allows the transmission of both cold and heat. The cold-hot pack remains flexible even when frozen and can be heated in the microwave.*

*Includes a fabric cover for applying it to the desired area.*

## COLD/HEAT PACK

## FEATURES

- **Non-toxic**
- **Includes protective fabric cover**
- **Individual packaging**
- **Microwaveable**
- **Flexible**

## SPECIFICATIONS

PRODUCT FAMILY: ORTOCELSOFT

PRODUCT NAME: Cold/Heat pack

INTENDED USE: To reduce pain and inflammation in post-surgical and post-traumatic treatments

Basic UDI-DI: 843603535SOFMP

### SIZE CHART

| Ref. | Measures (cm) |
|---|---|
| SOF0600 | 130 x 260 |

### COLOUR

- Blue

### COMPOSITION


**Bag:** Polyethylene/Polyamide
*LATEX FREE*

### COLD INSTRUCTIONS

- Place the gel pack into a freezer at least 2 hours before its use. Apply on the desired zone, wrapped in a dry towel or fabric cover. Keep the gel pack into the freezer to keep them cold for future applications.

### PACKAGING

- Case: 1 unit
- Box: 24 units

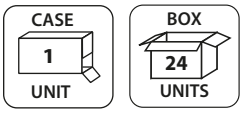


### HEAT INSTRUCTIONS

- Double boiler:
  Heat the gel pack in a container up to boiling point.
- Microwave: Heat the gel pack dipped in water during 5 minutes at maximum power of 600W.
  Apply on the desired zone, wrapped in a dry towel or fabric cover.

### ATTENTION

- Apply on the desired zone, wrapped in a dry towel or fabric cover. Is not recommended to use the gel pack directly on the skin
- Do not exceed more than 20 minutes in cold applications
- Patients with vascular insufficiencies must always consult a physician before using hot/cold therapy
- External use only
- Do not use a gel pack broken and throw it away
- Do not ingest
- Keep it away from children
- Non-toxic

MEDICARE SYSTEM
MEDI CARE SYSTEM S.L.U.
C/ Abat Escarré 15 - 08302 Mataró SPAIN
Tel.: +34 937 570 048 - Fax: +34 937 414 068
E-mail: info@medisystemsl.com · www.medicaresystem.es

CE Declaration according to Regulation (EU)745/2017 of the European Parliament & Council

*This data sheet could be changed without prior notice

CE