

SALVAFAST

ENGLISH

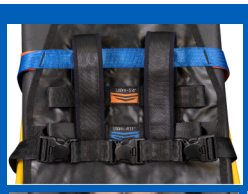

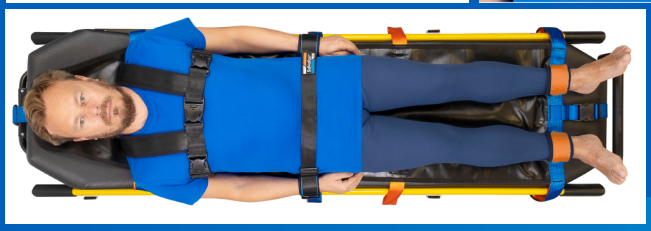


Salvafast Transfer *is an innovative device that provides quick restraint and safe transfer of altered, violent, aggressive or disoriented patients on emergency stretchers.*
It allows the transfer of patients from the emergency stretcher to another stretcher or bed, without the release of the restraint.
It consists of a resistant and durable mat provided of a set of harness, wristbands and ankle bands, with quick double Velcro® closure system and sliding base to improve the transfer.

SALVAFAST TRANSFER

FEATURES

One size fits all

Fits stretchers from 45 to 60 cm wide

Tested according to 5.4.1 of UNE-EN1865-1:2010+A1:2015
Patient handling equipment used in road ambulances - Part 1: General stretcher systems and patient handling equipment

Reinforced acetal double safety lock

Colour coded straps:
- Stretcher
- Handles
- Patient restraint

Foldable as roll for an easy storage and transport (Storage bag included)

SPECIFICATIONS

PRODUCT FAMILY: SALVAFAST

PRODUCT NAME: Salvafast Transfer

INTENDED USE: Fast restraint of disturbed, violent, aggressive or disoriented patients

Basic UDI-DI: 843603535SVTP8

SIZE CHART

Ref.	Measures (cm)		
	Mat	Wristbands	Ankle bands
SVT3710	172 x 45	15 - 30	20 - 40

COLOUR

Mat: Black, Orange

Straps: Blue, Orange, Black

COMPOSITION

Restraints: 80% PES - 20% Neoprene
Mat: 50% PVC - 50% PA
Side release buckles: Acetal (POM)
LATEX FREE

LAUNDRY INSTRUCTIONS

Wash the product by hand after each use with a damp cloth and a solution of ethanol or isopropanol. The product allows the use of surface disinfectants.

- Do not bleach
- Do not tumble dry
- Do not iron
- Do not dry clean

The product can be machine-washed up to 60ºC with a disinfectant detergent if necessary. However, it may worsen the product's deterioration.

DISINFECTION RECOMMENDATIONS

- Follow the disinfection guidelines of your health-care facility.
- Use a surface virucidal disinfectant in spray or wipes.
- You can use aqueous solutions of ethanol or isopropanol, quaternary ammonium (NH_4), or peracetic acid (CH_3CO_3H).
- You can use bleach occasionally. However, using sodium hypochlorite (HCL) solutions may reduce the product's lifespan.

PACKAGING

BULK PACKAGING

- 1 - Salvafast Transfer
- 1 - Storage bag

MEDICARE SYSTEM
MEDI CARE SYSTEM S.L.U.
C/ Abat Escarré 15 - 08302 Mataró SPAIN
Tel.: +34 937 570 048 - Fax: +34 937 414 068
E-mail: info@medisystemsl.com · www.medicaresystem.es

CE Declaration according to Regulation (EU)745/2017 of the European Parliament & Council
*This data sheet could be changed without prior notice

CE