MEDIROLLER

ENGLISH



Medi-roller Standard transfer board *allows the patient to be moved laterally comfortably and safely, without lifting them and using minimal force. Designed for transfers between beds, stretchers, operating tables, X-ray tables, etc.*

It is made up of a lightweight and robust polyethylene core covered with a polyester/PVC canvas sleeve welded by high frequency, achieving a compact, lightweight, and highly resistant surface. Tubular cover made of siliconized polyamide fabric, providing ultra-slippery properties.
The large foldable model has 2 handles for transport, and the short one-piece model has 1 handle.
Disposable cover available for infection prevention.

MEDI-ROLLER STANDARD

FEATURES



X-ray translucent and MRI safe



Ultra sliding sleeve cover



High frequency welded board



Lightweigth PE foam core



2 models:
Large - Foldable - 2 handles
Short - One piece - 1 handle

SPECIFICATIONS

PRODUCT FAMILY: Medi-roller

PRODUCT NAME: Medi-roller Standard

INTENDED USE: Perform transfers of patients in bed, stretcher, operating table or x-ray...

BASIC UDI-DI: 843603535TMRNE

SIZE CHART

Ref.	Model	Measures (cm)	Weight (Kg)
TMR6200	Foldable	180 x 50	2,00
TMR6220	Short	90 x 50	1,00

COLOUR

- Blue board
- Orange cover

COMPOSITION

Board: Ext: 70% PVC - 30% Polyester
Int: Polyethylene
Cover: Silicone coated Polyamide
LATEX FREE

LAUNDRY INSTRUCTIONS

Board: Clean with damp cloth and mild soap or surface disinfectant detergent.

Cover: - Machine wash max. 60ºC - Do not tumble dry
- Do not bleach - Do not iron
- Do not dry clean




DISINFECTION RECOMMENDATIONS

- For infection prevention and control use the disposable protective cover (ref. TMR6270).
- Use a surface disinfectant or disinfectant wipes if required. For cover disinfection machine wash at 60ºC using a disinfectant detergent. Do not use softener, chlorine nor whitening concentrates.

WARNING

- Maximum weight: 170 Kg.
- Maximum distance between surfaces: 20 cm

PACKAGING

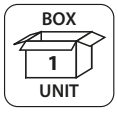

- 1 - Medi-roller Standard board
- 1 - Medi-roller Standard cover
- 1 - Medi-roller Standard disposable cover

MEDICARE SYSTEM
MEDI CARE SYSTEM S.L.U.
C/ Abat Escarré 15 - 08302 Mataró SPAIN
Tel.: +34 937 570 048 - Fax: +34 937 414 068
E-mail: info@medisystemsl.com · www.medicaresystem.es

CE Declaration according to Regulation (EU)745/2017 of the European Parliament & Council

*This data sheet could be changed without prior notice

CE